




Start Sensor Session



These instructions are specific to Dexcom G7.

For information on starting a Dexcom G6 sensor session, please refer to the user guide.*

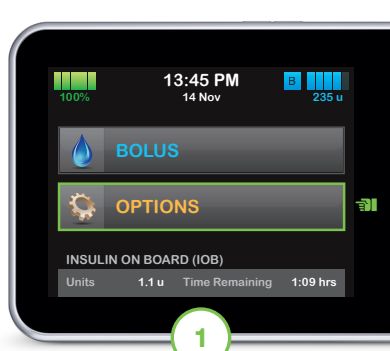

From the Home screen, tap **OPTIONS**.

Note: These instructions are provided as a reference tool for pump users and caregivers who are already familiar with the use of an insulin pump and with insulin therapy in general. Not all screens are shown. For more detailed information on the operation of the t:slim X2 pump, please refer to its user guide.



Tap the down arrow and then tap **My CGM**.



If this is the first time using a Dexcom G7 sensor, tap **Select Sensor** and then tap **Dexcom G7**. Otherwise, skip to Step 4.




Pump Tip: The Dexcom G7 sensor needs a 30-minute startup period, which begins immediately after insertion, before trend information displays on the pump screen. When pairing your G7 sensor, make sure all other sensors are at least 6 meters (20 feet) away from your pump.



Tap **START G7 SENSOR**.


Note: If the user utilizes the Dexcom G7 mobile app, the user must start a sensor session and enter the pairing code on both the Dexcom G7 mobile app and t:slim X2 insulin pump.

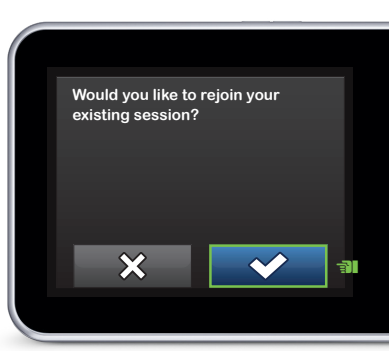


If the user is continuing with an active sensor session and is not using the Dexcom G7 mobile app, tap  and then skip to Step 7. Otherwise, go to Step 5 to continue.



Enter the code found on the side of the sensor applicator and tap  to continue.

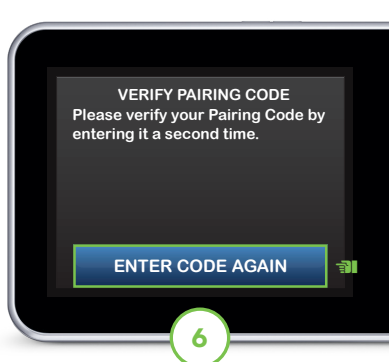


Tap **ENTER CODE AGAIN** to verify the pairing code and then tap  to continue.

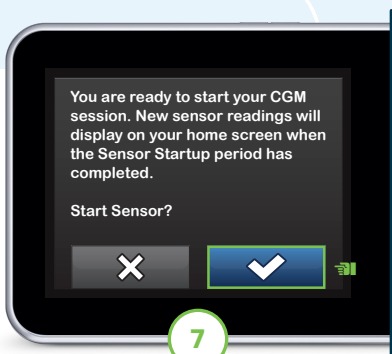


Tap  to start the sensor session startup period.

The SENSOR STARTED screen will appear to confirm.

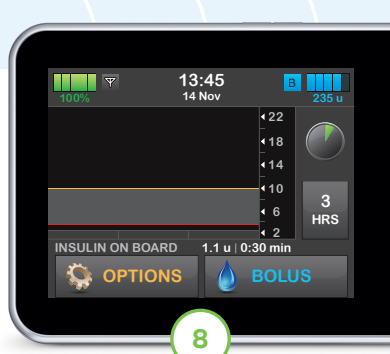

The pump will return to the Home screen with a CGM trend graph and sensor startup countdown symbol.

Note: If continuing an active sensor session, the startup period will take five to 10 minutes. For a new session, it may take up to 30 minutes.

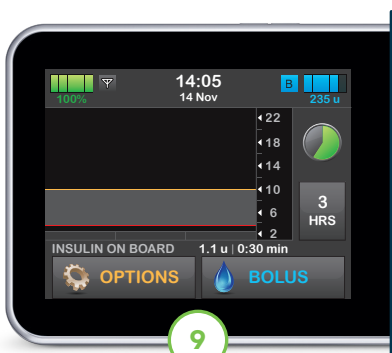

The symbol fills in over time to show how much time is left.

Note: During this startup period, the pump will not display sensor data or allow use of predictive technology. Users should continue to use a blood glucose (BG) meter and test strips in order to make treatment decisions.




After the startup period, the new active sensor session will begin.

Note: If you see a Sensor Not Started alert, contact Dexcom support. If glucose alerts and readings do not match symptoms or expectations, use a BG meter to make diabetes treatment decisions.

Connection Tips

If the Out of Range icon is on the insulin pump screen (pictured below), then the sensor is not communicating with the pump. When a sensor glucose reading cannot be provided, “- -” shows in the place where the sensor glucose value is normally shown.

- ✓ For ideal connectivity, it is suggested that the user faces the pump screen out and away from their body and wear the pump on the same side of the body that they wear the sensor
- ✓ Make sure that nothing is rubbing or obstructing the sensor
- ✓ Please note that wireless communication does not work well through water so the transmission range is lower

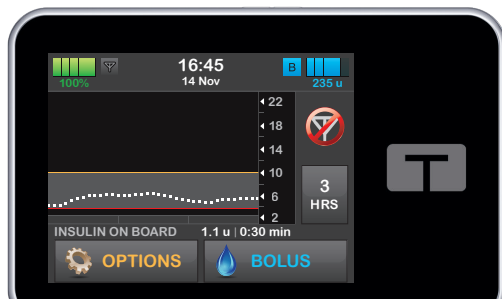

* CGM sold separately.

Important Safety Information: The t:slim X2 insulin pump is intended for the subcutaneous delivery of insulin, at set and variable rates, for the management of diabetes mellitus in persons requiring insulin. The pump is able to reliably and securely communicate with compatible, digitally connected devices. The pump is indicated for use in individuals 6 years of age and greater. The pump is intended for single user, home use and requires a prescription in certain countries. The pump is indicated for use with NovoRapid or Humalog U-100 insulin. The System is not indicated for use in pregnant women, people on dialysis, or critically ill users. Users of the pump and the System must: be willing and able to use the insulin pump, CGM, and all other system components in accordance with their respective instructions for use; test blood glucose levels as recommended by their healthcare provider; demonstrate adequate carb-counting skills; maintain sufficient diabetes self-care skills; see healthcare provider(s) regularly; and have adequate vision and/or hearing to recognize all functions of the pump, including alerts. The t:slim X2 pump must be removed before MRI, CT, or diathermy treatment. Visit tandemdiabetes.com/safetyinfo for additional important safety information.

© 2024 Tandem Diabetes Care, Inc. All rights reserved. Tandem Diabetes Care, Tandem logos, and t:slim are registered trademarks of Tandem Diabetes Care, Inc. in the United States and/or other countries. Dexcom, Dexcom G6, Dexcom G7, and any related logos and design marks are either registered trademarks or trademarks of Dexcom, Inc. in the United States and/or other countries. All third-party marks are the property of their respective owners. ML-1011948_B

